




From modelling to McStas and more

A short guide

E. Farhi, *ILL*

Data formats

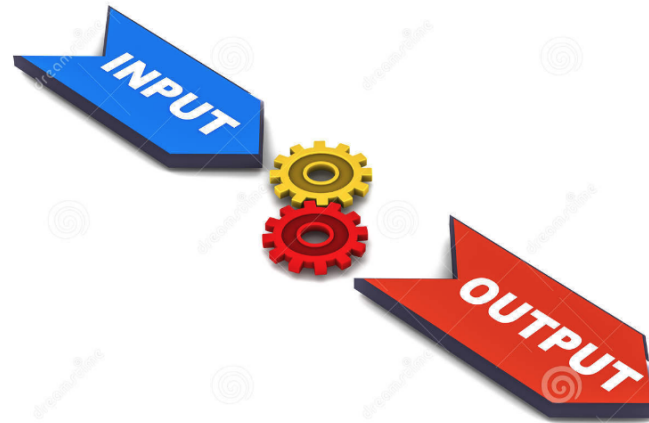

McStas Data formats



McStas uses a number of input data formats, all text based.

- **Lau** (Single_crystal, PowderN, Isotropic_Sqw)
- **Laz** (PowderN, Isotropic_Sqw)
- **Sqw** (Isotropic_Sqw)
- **qSq** (Isotropic_Sqw) – mostly for liquids
- **nxs** (Sample_nxs) – this is ***not*** a NeXus/HDF
- **PDB** (SANSPDB)




The question is how to generate these files...

By hand ?




Material modelling formats

- MD (classical or DFT)
 - trajectories, many data formats
- Lattice dynamics:
 - PhonoPy FORCE_SETS
 - ...

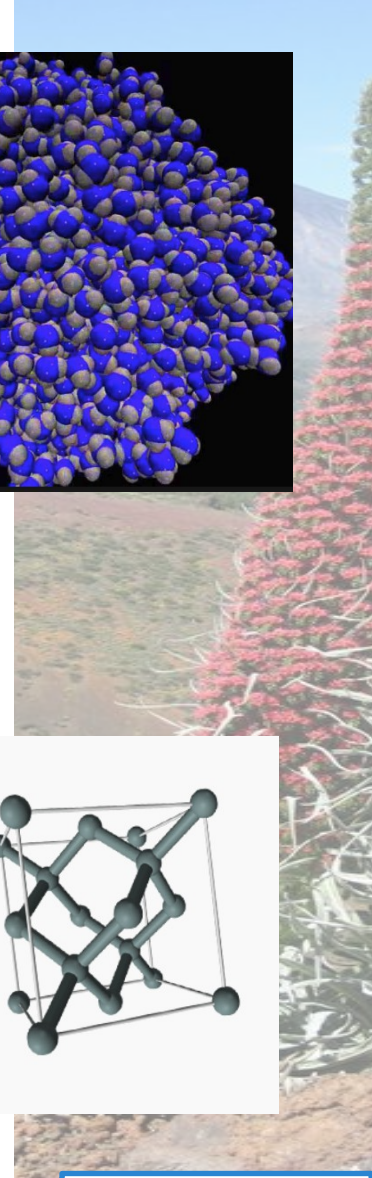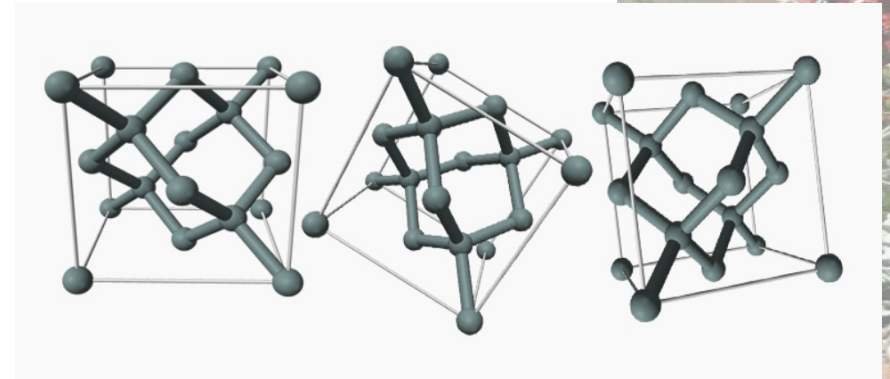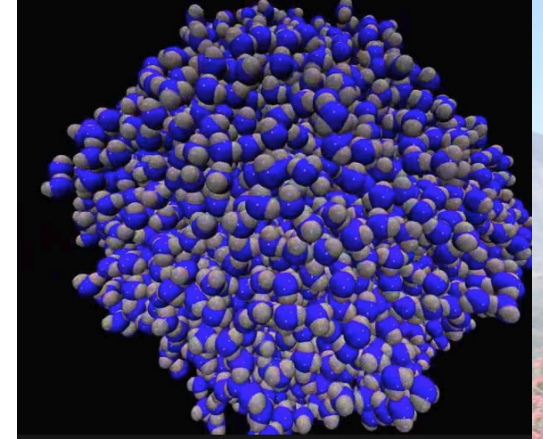

Nuclear data base formats


- ENDF, JEFF, JENDL, ...
 - all use a standard format for a set of materials. Most are monoatomic / incoherent
- ACE
 - Same as ENDF, but in binary for MCNP

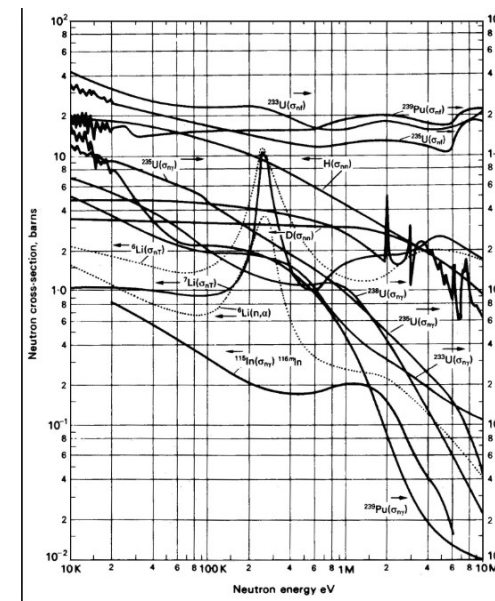



Experiments

- Intensity(d) \rightarrow $F(Q)$ \rightarrow $S(Q)$ and Laz
- Intensity(HKL) \rightarrow $F_2(\text{HKL})$ \rightarrow Lau
- Intensity(theta,ToF) \rightarrow $S(q,w)$
- Often saved as text and/or NeXus

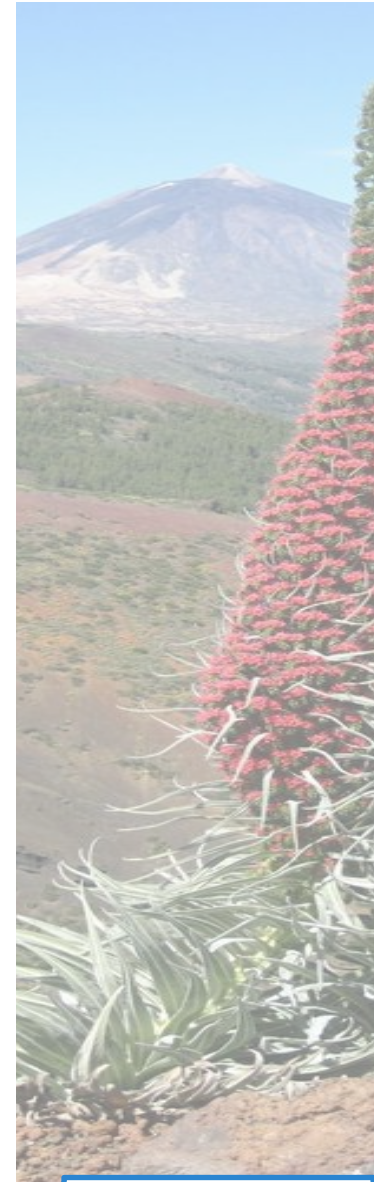

Tools

- iFit
- Mantid
- Python/ASE (IO)
- MDANSE/nMoldyn (converters)
- Scripts of all kinds
 - OpenBABEL, ...
- cif2hkl
- bash






Material DataBases


- ICSD <<http://icsd.ill.fr/icsd/index.html>>
- COD <<http://crystallography.net>>
- <<https://www.materialsproject.org/>>
- PhononDB <<http://phonondb.mtl.kyoto-u.ac.jp/>>
- Nomad <<https://repository.nomad-coe.eu/>>

Mostly provide structures...

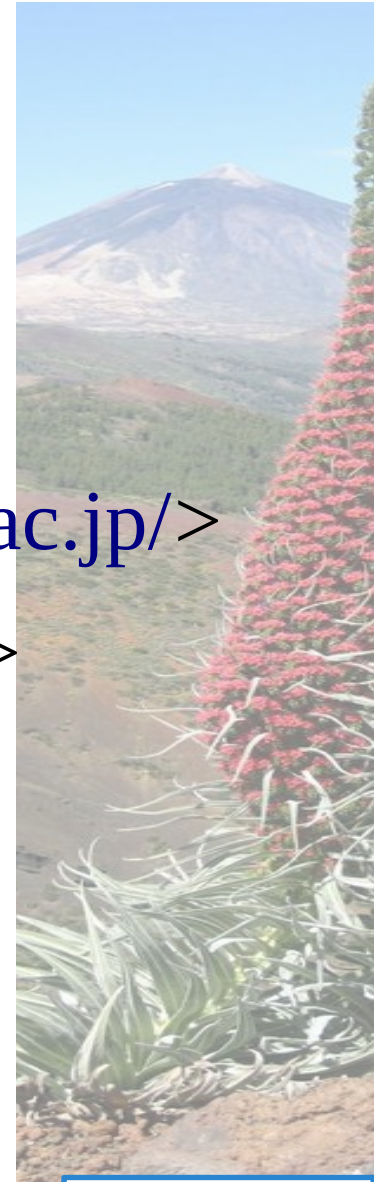

McStas: structure: powder

- For powders, the pure elastic diffraction is easy to prepare: LAZ
 - PowderN
 - Single_crystal(powder=1, ...)
 - Isotropic_Sqw

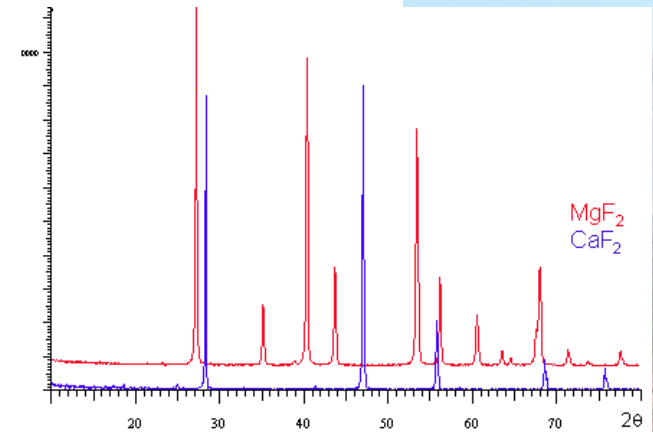

CIF file
FullProf/ShelX


Terminal:
`cif2hkl --powder file.cif`

McStas Laz file

CIF file
FullProf/ShelX

iFit:
`cif2hkl('file.cif', 'p')`

McStas Laz file

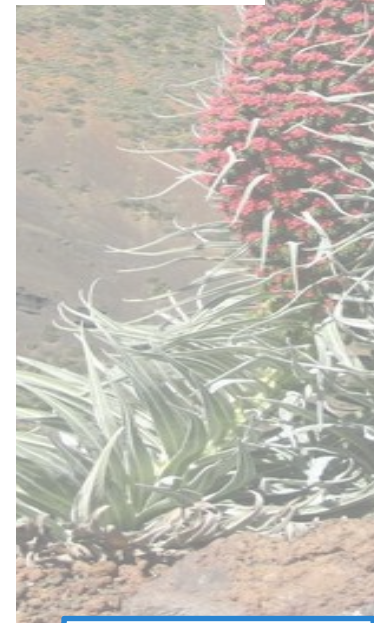



McStas: structure: SX

- For single crystals, the pure elastic diffraction is easy to prepare: LAU
 - PowderN
 - Single_crystal
 - Isotropic_Sqw

CIF file
FullProf/ShelX

Terminal:
`cif2hkl --xtal file.cif`

McStas Lau file

CIF file
FullProf/ShelX



iFit:
`cif2hkl('file.cif', '', 'x')`

McStas Lau file



McStas: structure: COD


- Cif2hkl from iFit can input chemical formula
- Requires internet connection to COD

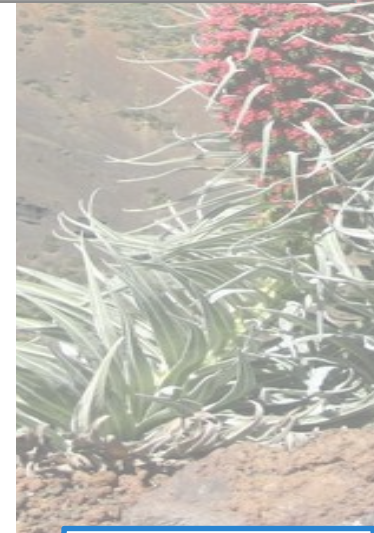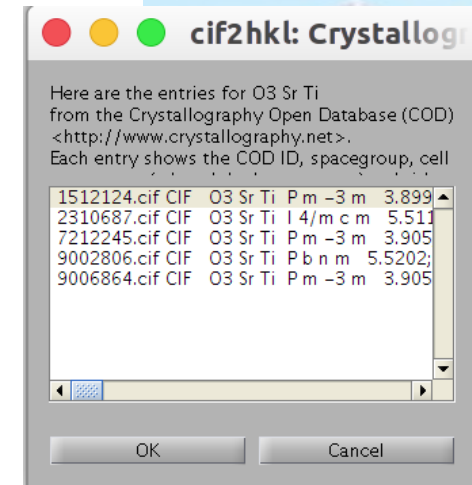



McStas: structure: others

- Other McStas structure file formats
 - SANSPDB and similars: get from Data Bases
 - Sample_nxs: manual edit only (nxs, text)
 - Isotropic_Sqw: manual edit only (qSq, text)
 - NCrystal_sample*: (ncmat,nxs,lau)

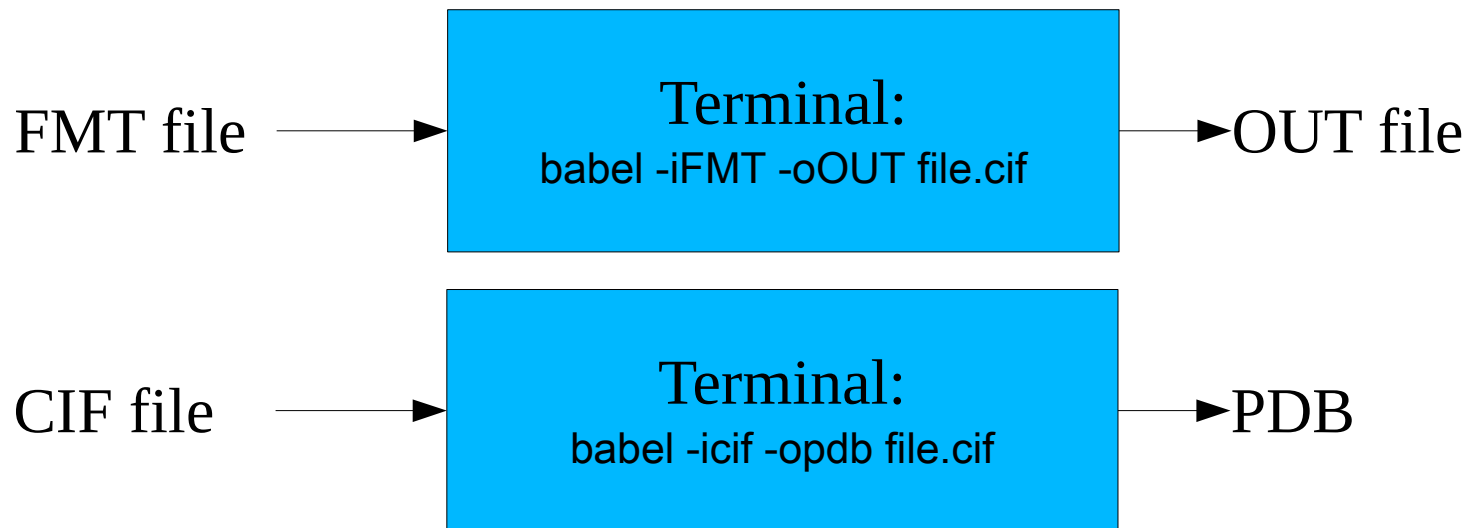

OpenBabel supports many data files

*in McStas from v. 2.5





McStas: dyn from ToF spec

- Import and correct data with e.g. LAMP, Mantid
- Integrate radially ($|q|$) and compute $S(q, \omega)$
- Export to e.g. HDF/NeXus or other
- Import then export with iFit:
 - `sqw=iData_Sqw2D('file')`
 - `sqw.saveas('file.sqw', 'mcstas')`
- Can/should always edit manually (text)

```
# weight    4.002  in [g/mol]
# density   0.4784 in [g/cm^3]
# sigma_abs 0.00747 absorption scattering cross section in [barn]
# sigma_coh 1.34   coherent scattering cross section in [barn]
# sigma_inc 0      incoherent scattering cross section in [barn]
# Temperature 2      in [K]
# classical 0      experimental, contains Bose factor
```





McStas: dyn from exp. DOS

- Import and correct data with e.g. LAMP, Mantid
- Estimate $g(w)$ with e.g. LAMP, Mantid, MUPHOCOR
- Export to e.g. HDF/NeXus or other
- Import then export with iFit:
 - `g=iData_vDOS('file')`
 - `inc=g.incoherent;`
 - `inc.saveas('inc.sqw','mcstas')`
- Can/should always edit manually (text)





McStas: dyn from nuclear DB

- Import and export with iFit:
 - `sqw=iData_Sqw2D('tsl.endf')`
 - `saveas(sqw(1), 'inc.sqw', 'mcstas')`
- Can/should always edit manually (text)

```
# weight    4.002  in [g/mol]
# density   0.4784 in [g/cm^3]
# sigma_abs 0.00747 absorption scattering cross section in [barn]
# sigma_coh 1.34   coherent scattering cross section in [barn]
# sigma_inc 0      incoherent scattering cross section in [barn]
# Temperature 2      in [K]
# classical 0      experimental, contains Bose factor
```





McStas: dyn from MD

- Import trajectory into MDANSE (converters)
- Double click the trajectory to allow analysis
- Use Plugins:
 - Analysis/Scattering/DCSF and DISF
 - *Analysis/Scattering/Structure Factor*
 - *Analysis/Dynamics/Density of States*
- Import then export with iFit:
 - `sqw=iData_Sqw2D('DCSF.nc')`
 - `sqw.saveas('coh.sqw','mcstas')`
- Can/should always edit manually (text)





McStas: dyn from nuclear DB

- Import and export with iFit:
 - `sqw=iData_Sqw2D('tsl.endf')`
 - `saveas(sqw(1), 'inc.sqw', 'mcstas')`
- Can/should always edit manually (text)

```
# weight    4.002  in [g/mol]
# density   0.4784 in [g/cm^3]
# sigma_abs 0.00747 absorption scattering cross section in [barn]
# sigma_coh 1.34   coherent scattering cross section in [barn]
# sigma_inc 0      incoherent scattering cross section in [barn]
# Temperature 2      in [K]
# classical 0      experimental, contains Bose factor
```





McStas: dyn from Lattice Dynamics



- Using PhonoPy:
 - iFit: `s4d=sqw_phonons('phonopy_directory')`
- Computing on the fly:
 - iFit: `s4d=sqw_phonons('cod: La Mn O3')`
 - `p=powder(s4d)`
 - `d=iData_Sqw2D(iData(p, [], 0:10, 0:100))`
 - `d.saveas('coh.sqw', 'mcstas')`
- Other:
 - Import in iFit, convert to Sqw2D, export.





Neutron events

- McStas handle ‘rays’ as packets of neutron events.
- Summing up builds intensity and distributions.
- McStas can read/write ‘events’ from other software:
 - MCNP TRAC and SSW
 - TRIPOLI
 - MCPL
- McStas can generate events as well for Mantid via NeXus/HDF files.

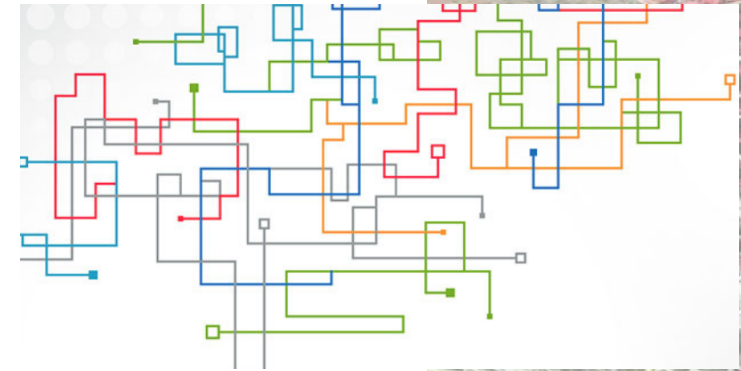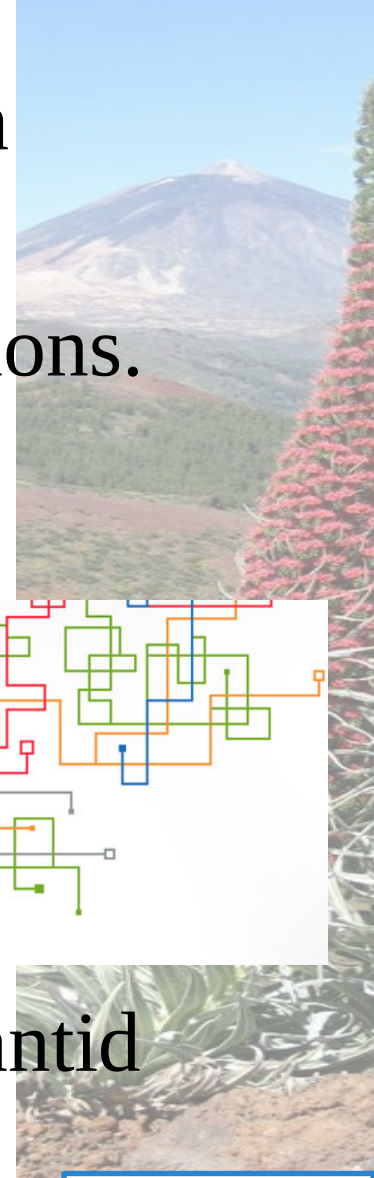

McStas automation: bash



- Can use cif2hkl and mcrun commands:
 - **cif2hkl** 'file' --out 'file.laz' --powder
 - **mcrun** --dir=test **instr** Powder=file.laz
 - **mcplot.pl** -psc test

```
[~]$ bash
```



McStas automation: iFit



- iFit is a generic infrastructure which gathers data sets handling, fitting, $S(q,w)$ models, and McStas hooks.

iFit





McStas iFit: build

- McStas can be controlled from within iFit.
- Open Matlab/iFit
- Create the McStas model with:
 - `model = mccode('instr')`
- Plot the geometry with:
 - `plot(model)` *% has contextual menus*
- Edit the instrument and re-compile
 - `edit(model)`

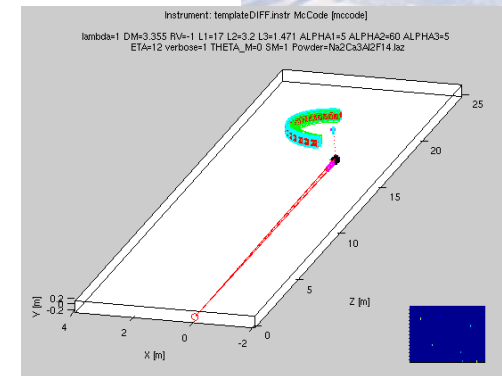



Figure 2: TextEdit: /tmp/tp109e7fc7_fd80_

```

File Edit Help
McStas instrument definition URL=http://www
*****
* Instrument: Template monochromator Diffractometer
*
* %Identification
* Written by: E. Farhi
* Date: 13 Apr 2006
* Origin: LLB/ILL
* Release: McStas CVS_080624
* Version: $Revision$
* %INSTRUMENT_SITE: Templates
*
* Simple monochromator Diffractometer for powders
*
* %Description
  
```

McStas iFit: eval


- Run with (default 1e6 event and pars)
 - `data = iData(model, [], nan);`
 - `subplot(model) % plot monitors, contextual menus`
- Specify parameters
 - `data = iData(model, 'lambda=2.36; coh=Cu.laz')`
- Do a scan:
 - `data = iData(model, 'lambda=[1.2 2.4 3.6]')`
- Change neutron events #
 - `model.UserData.options.ncount = 1e7;`

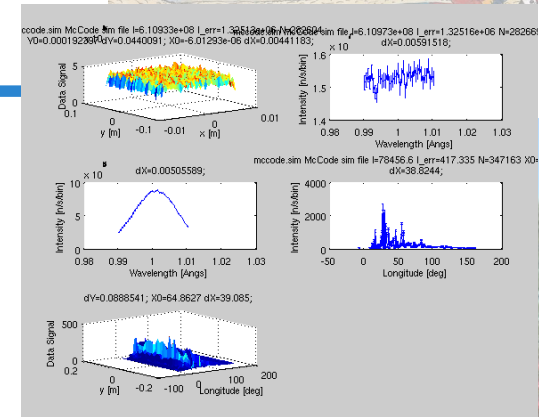



McStas iFit: optim

- Fix all parameters but *lambda*, Maximize model value:
 - `mlock(model, 'all');` `munlock(model, 'lambda')`
 - `xlim(model, 'lambda', [1 3]);` % bounds
 - `fmax(model, 'lambda=2.36', '', nan)`
 - `fmax(model, 'lambda=2.36', 'OutputFcn=fminplot', nan)`

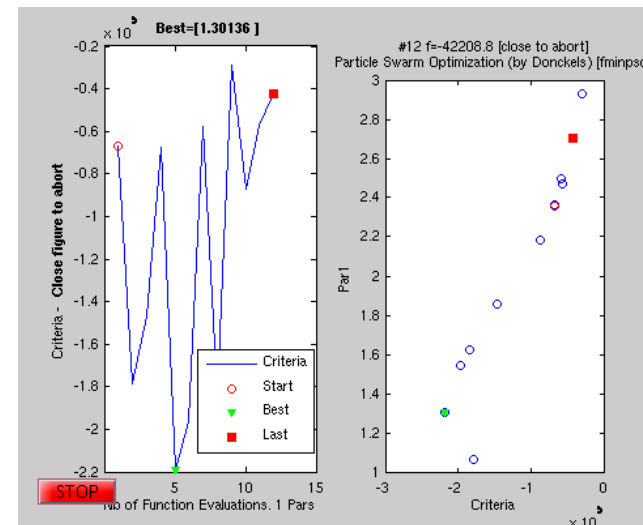



McStas iFit: advanced

- You can add McStas models:
 - `model=mccode('instr1')+mccode('instr2') + ...`
- For instance:
 - Instr1: structure
 - Instr2: spin-wave
 - Instr3: phonons
 - Instr4: incoherent
- The different models can also be assembled as a set of samples chosen at execution in a single McStas instrument.

